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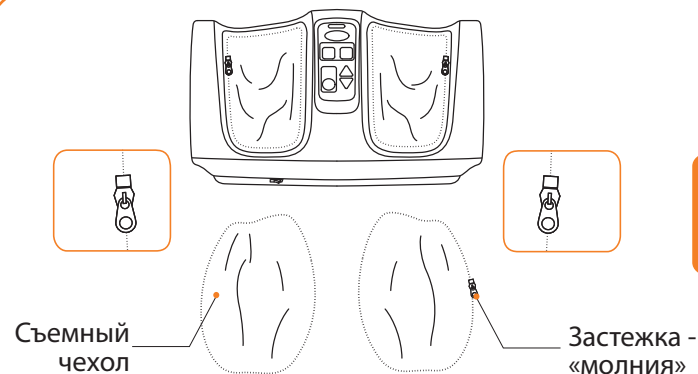
## CONTENTS

Features.....	15
Important safety instructions.....	16
Operating instruction .....	18
Working modes.....	20
Troubleshooting.....	22
Care & maintenance .....	23
Specifications.....	23

## FEATURES

- Deep kneading massage for feet and calves
- For home and office
- Rollers' adjustable speed and rotation direction
- 2 automatic massage modes
- Timer: 15 minutes
- Two extra pairs of easy-to-clean covers included

## IMPORTANT SAFETY INSTRUCTIONS

**Please read this user's manual attentively.**

- Do not put hands or fingers near the massage mechanism or rollers while this product is running.
- Touching the mechanism can result an injury from the rollers squeezing together.
- If the supply cord is damaged, it must be replaced by the manufacturer or its service agent or a similarly qualified person in order to avoid a hazard.
- The appliance should never be left unattended when plugged in. Unplug from the outlet when not in use, and before putting on or taking off parts.
- If you feel any abnormal discomfort during the operation of this product, turn off the power immediately and do not use the product. Consult your physician.
- Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return this product to a service center for examination and repair.
- Always keep the power cord away from heated surfaces.
- Never drop or insert any object into any opening.
- When using this product after it has been in storage, check to see that the product is operating properly prior to use.
- To disconnect, turn all controls to the OFF position, then remove plug from the outlet.
- Do not massage feet or calves that are swollen, inflamed or covered with skin eruptions.
- Connect this appliance to a properly grounded outlet only.
- Do not use this product if noise above normal sound is heard.
- Do not stand or put heavy objects on the massager. Use only while seated.
- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety. Children should be supervised to ensure that they do not play with the appliance.

**SAVE THESE INSTRUCTIONS**

**ATTENTION!**

- **Do not** operate where aerosol (spray) products are being used or where oxygen is being administered.
- **Do not** operate under blanket or pillow.
- **Do not** stand on the massager.
- **Do not** pull by the AC cord.
- **Do not** use outdoors.
- **Do not** carry this appliance by the power cord or use the cord as a handle.
- **Do not** reach for a product that has fallen into water. Unplug immediately.
- **Do not** use in the place with high humidity.

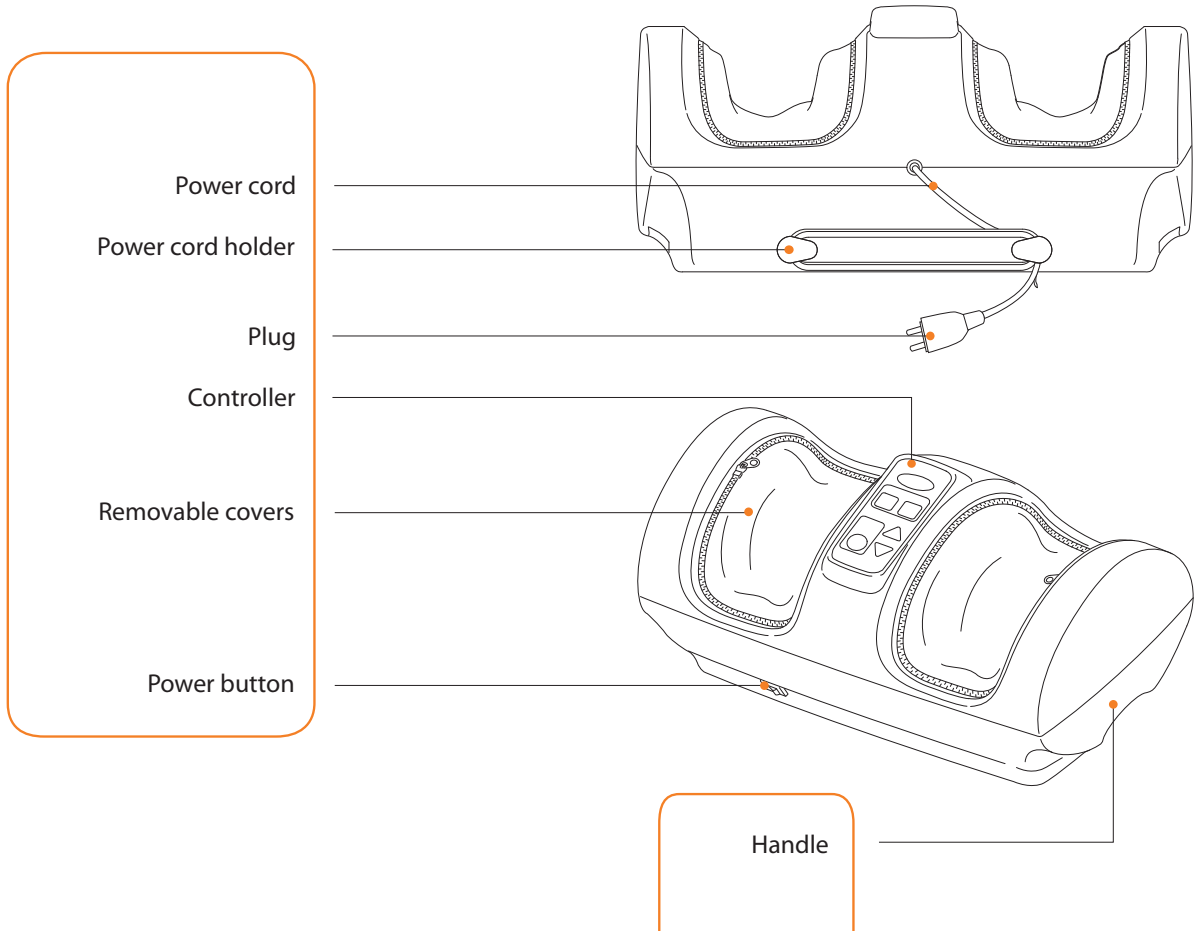
**Persons who are undergoing medical treatment as well as the following persons should consult a physician before use:**

- Persons who have injured their legs or feet in the past, persons who have numbness in their feet.
- Persons having a bone condition, such as osteoporosis.
- Persons having a circulatory disorder (heart disease, angiopathy, hypertension, etc.).
- Persons having a respiratory organ disorder.
- Persons using a body implanted medical electronic device.
- Persons having malignant tumors.
- Persons having symptoms of a circulatory disorder (thrombosis, serious aneurysm, chronic varicosis, etc.), dermatitis, or infectious skin diseases.
- Persons having a sensory disorder due to an advanced peripheral circulatory disorder, such as that resulting from diabetes.
- Persons having skin injuries.
- Persons having a fever (38°C or higher) due to disease.
- Persons who are pregnant or thought to be pregnant or who are menstruating.
- Persons who feel physical discomfort and require rest.
- Persons wishing to use the unit for rehabilitation purposes.
- Persons other than the above who do not feel physically normal.


Consult your physician  
before use!


# OPERATING INSTRUCTION

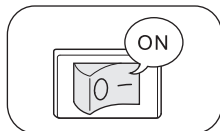- Place the massager on a flat surface.
- Plug the power cord into the 220 V socket.
- Press the power button on the front part of the massager.

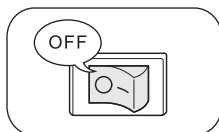- Place your feet into the massager.
- Press **On/Off** button on the controller. The massage rollers will start movement.
- After 15 minutes of the massage, the massager will turn off automatically.
- After finishing the massage, press the power button on the front part of the massager.


- Unplug the massager.

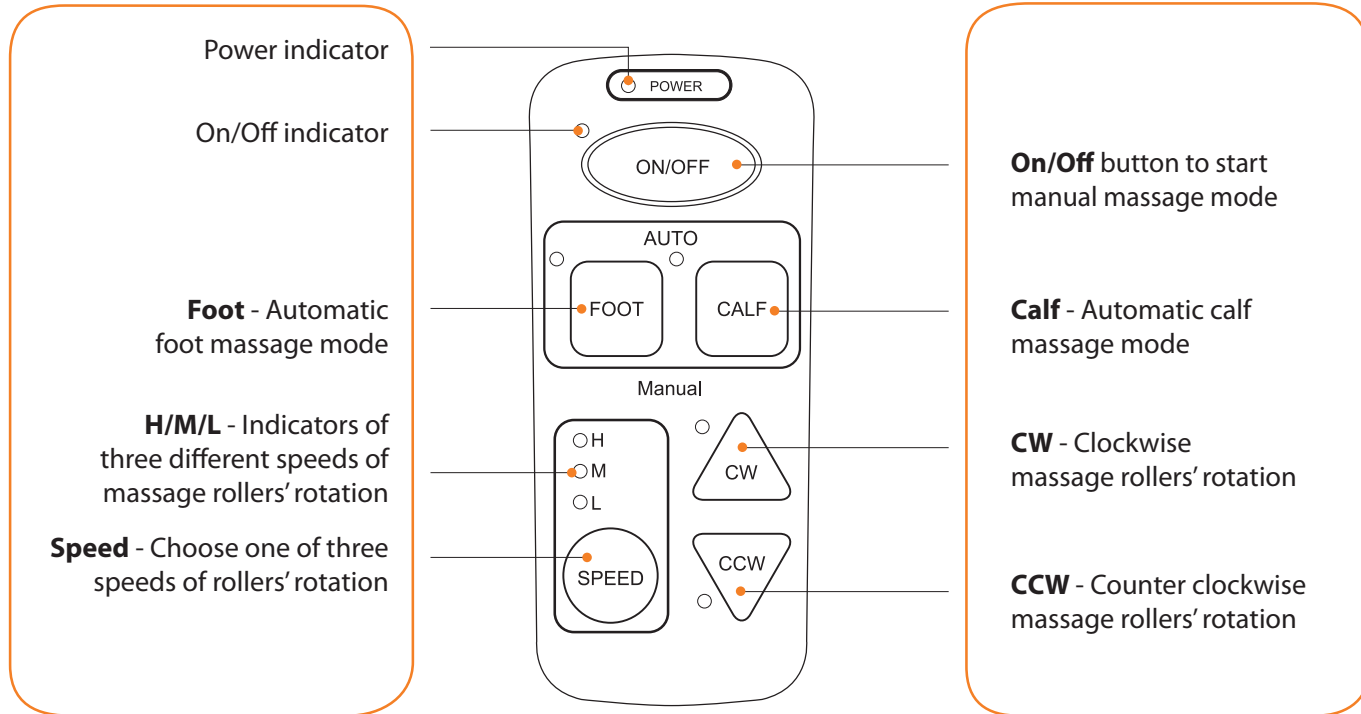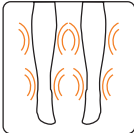
In case you need to stop the massage immediately, press **On/Off** button on the controller.

The massager working time is 15 minutes at one session.

# WORKING MODES

## Controller

MODE	DESCRIPTION	BENEFITS
<p><b>FOOT MASSAGE</b></p> 	<p>Press <b>Foot</b> button on the controller to start an automatic foot massage mode. After you press this button the corresponding indicator will light up.</p>	<p>Relieves tension and soreness in your feet by lifting and stretching muscles.</p>
<p><b>CALF MASSAGE</b></p> 	<p>Press <b>Calf</b> button on the controller to start an automatic calf massage mode. After you press this button the corresponding indicator will light up.</p>	<p>Improves circulation and relieves muscle tension.</p>
<p><b>MANUAL MODE</b></p>	<p>After you turn the massager on by pressing the power button, press <b>On/Off</b> button on the controller to start the manual massage mode. During this mode you can choose one of three different speeds of massage rollers' rotation: high/medium/low by pressing repeatedly <b>Speed</b> button on the controller.</p> <p>During the manual massage mode you can also choose the direction of massage rollers' rotation: clockwise or counter clockwise.</p> <p>Press one of the buttons <b>CW</b> or <b>CCW</b> on the controller:</p> <p><b>CW</b> — clockwise.</p> <p><b>CCW</b> — counter clockwise.</p>	<p>Stimulates your entire body through the numerous acupressure points located on the surface of your feet. Foot massage calms down your nervous system and improves the blood circulation.</p>

## TROUBLESHOOTING

PROBLEM	CAUSE	SOLUTION
Massage stopped	The preset maximum period of operation has elapsed	Let the device cool down for at least 15 min before using again. Switch the device on again
	The connection to the mains has been interrupted	Check the plug connection and the contact of the main plug with the main supply
	<b>On/Off</b> button on the controller is not pressed	Press <b>On/Off</b> button on the controller
	Power button is not pressed on the front part of the massager	Press power button on the front part of the massager
Loud noise occurs during the massage	The noise can be caused by the massage mechanisms' movement	This situation is normal and it is not a trouble of the massager

## CARE & MAINTENANCE

### Storage

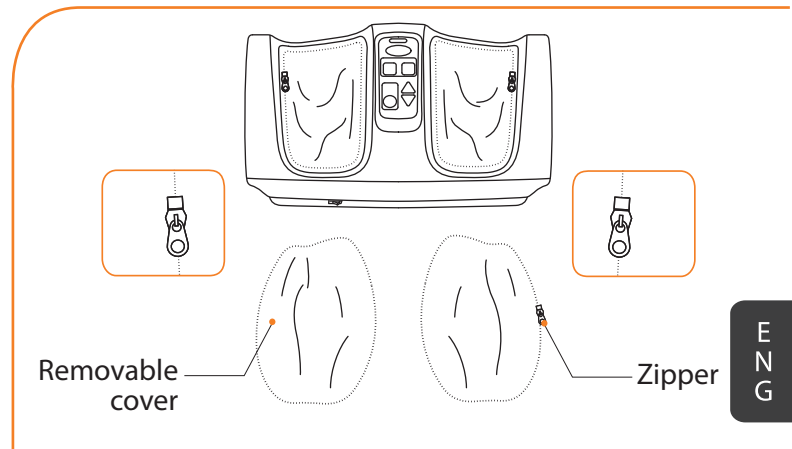
- Do not wrap the power cord around the device before storing.
- Store the foot massager in its original box and in a safe, dry and cool place.
- Store it out of the reach of children.
- Avoid contact with sharp edges or pointed objects as this might damage the fabric of the covers.

### Cleaning

- Unplug the device.
- Use a soft, damp cloth to clean the device.
- To clean the removable covers of the massager, unzip them (see the picture to the right), take them off and wash in a delicate mode.
- Use a soft brush to clean dirty connections and contacts.
- Dry the device with a soft cloth to minimise potential water damage.
- Wait until the device is completely dry before using it again.

## SPECIFICATIONS

Product name	Foot massager
Model	ANGEL FEET
Voltage	220 V, 50/60 Hz
Power	45 W
Timer	15 minutes

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US Medica LLC., 1201 N.Orange St, Wilmington, DE 19801, USA